

iBrin-3000E

Electronic Touch Screen Brinell Hardness Tester



Video



Contact us

Mikrosize Precision Instrument Co.,Ltd

A-4035 RuiFeng Business Expo, Wuhu City, China , 241000.

Web: www.mikrosize.com

Email: mikrosize@mikrosize.com



Feature and Application

Product Feature

- Closed loop loading with high-precision load cell and single - chip microcomputer control ensure precise automatic loading and unloading, guaranteeing reliable test results.
- Variety of test force levels, ranging from 0.613 KN to 29.42 KN, meeting the hardness testing requirements for different materials.
2.5mm,5mm,10mm three carbide brinell ball indenters
- The test force dwelling time can be flexibly adjusted from 1 to 99 seconds, adapting to the characteristics of various materials and testing standards.
- The maximum test height is 280mm, and test throat is 150mm, providing ample operating space.
- Instrument function on touch screen is clearly divided, operating simple: data viewing, processing, and calibration can be easily completed.
- iBrin-3000E complies with standards such as GB/T231.2 - 2002 ensuring accuracy.

Product Application

- iBrin-3000E brinell hardness tester is a bench - type machine, which is used to test the brinell hardness of annealed and normalized workpieces, cast parts, non - ferrous metals, various relatively soft components, or unquenched steel parts.
- It is suitable for use in metrology, metal metallurgy, chemical industry, machinery manufacturing, industry, as well as scientific research institutions of universities and colleges.
- Accurately detect the hardness of annealed and normalized parts to control the quality of metal processing.
Effectively measure the hardness of cast parts to ensure stable performance of products.
- Reliably test the hardness of non - ferrous metals to support the R & D and innovation of materials.
Precisely measure the hardness of soft components to ensure a good fit during assembly.
- Quickly identify the hardness of unquenched steel parts to improve production efficiency.




Product Details



Product Details



Touch Screen Function

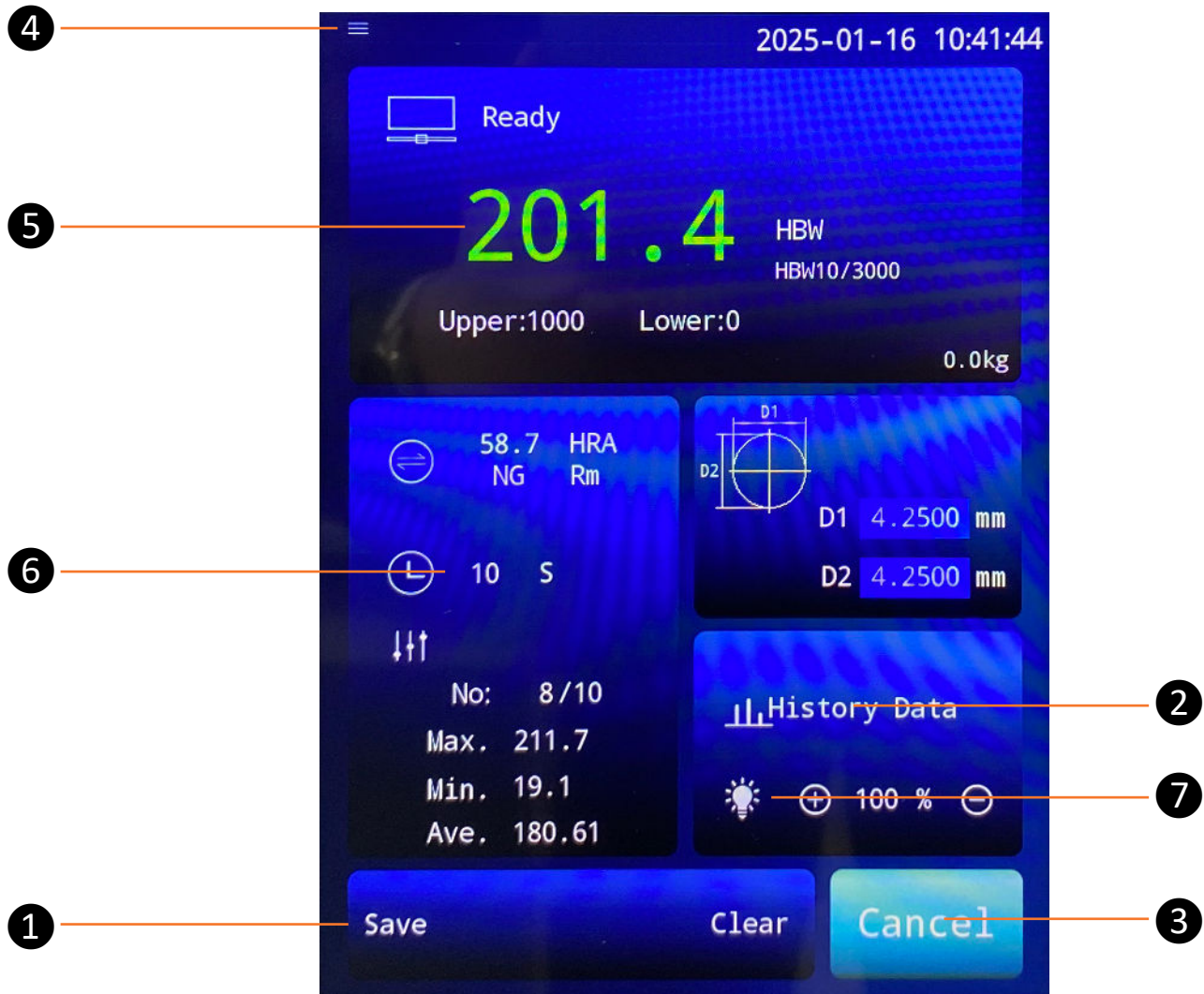

- 1.Save: Save current test data and store it in the historical data.
- 2.History Data: View the previously saved historical data. You can scroll up and down to check, and also delete, print or export the historical data to a USB.
- 3.Cancel: Press "Cancel" during operation to stop the running process, and the indenter will return to initial state.
- 4.Function keys: Used for force value calibration, encoder calibration, and language switching.
- 5. Hardness Test Result
- 6. Dwelling Time Setting
- 7. Illumination Adjusting


Product Details

Brinell Hardness Scale Selection

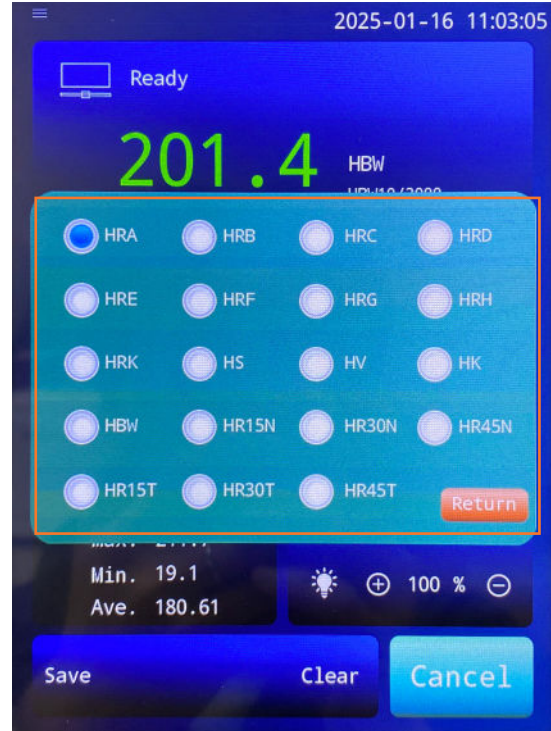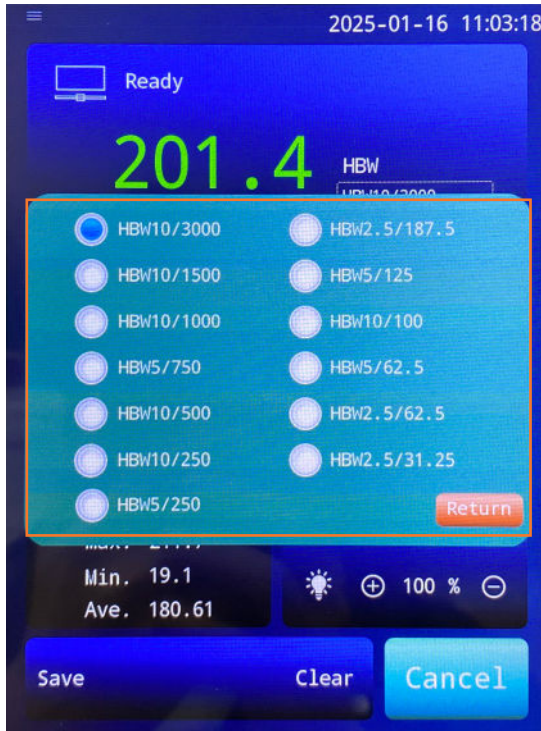

Conversion Scale Selection

Measured Data Input

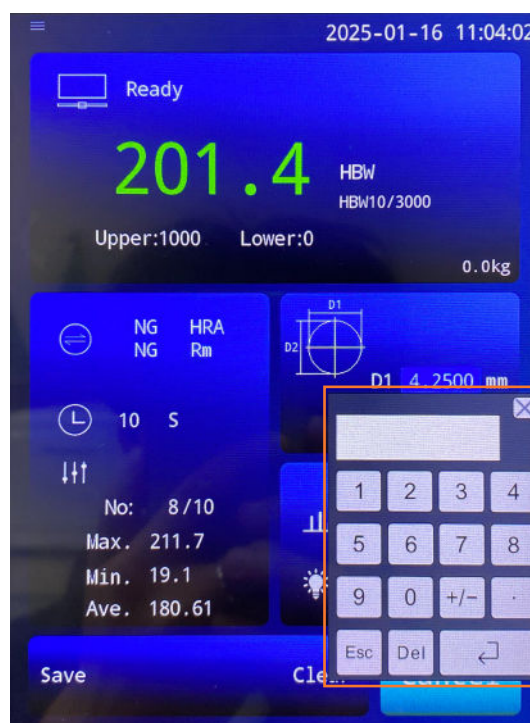

Product Details

BRM-2 Portable Brinell Reading Microscope



- BRM-2 brinell reading microscope with light source is specially used for measuring the indentation size of Brinell hardness. Its own light source ensures that the indentation image and scale line are super clear in the environment with no light, makes indentation measurement is easy and accurate.

Model	BRM-2
Total Magnification	20X
Eyepiece Microscale Range	0 mm to 6 mm
Eyepiece Maximum field of View	9 mm
Effective Measuring Range	0 mm to 6 mm
The Minimum Reading of Microdrum	0.01mm
Measurement Accuracy	+ / - 0.01mm
Weight	500g
Dimensions	82X35X155mm

Technical Specification

Product Name		Electronic Touch Screen Brinell Hardness Tester
Model		iBrin-3000E
Data Display&Output		LCD touch screen
Test Force	kg	62.5kgf、 100kgf、 125kgf、 187.5kgf、 250kgf、 500kgf、 750kgf、 1000kgf、 1500kgf、 3000kgf
	N	612.9N, 980.7N, 1226N, 1839N, 2452N, 4903N, 7355N, 9807N, 14710N, 29420N
Test Force Resolution		0.1Kg
Brinell Scale		HBW2.5/31.25、 HBW2.5/62.5、 HBW2.5/187.5、 HBW5/62.5、 HBW5/125、 HBW5/250、 HBW5/750、 HBW10/100、 HBW10/250、 HBW10/500、 HBW10/1000、 HBW10/1500、 HBW10/3000
Test Force Accuracy		62.5 ~ 250Kgf≤1%, 500 ~ 3000Kgf≤0.5%
Reading Microscope		Resolution :0.01 mm (The minimum reading of microdrum)
Ball Diameter		Φ2.5mm、 Φ5mm、 Φ10mm
Test Space		Test Height: 280mm Test Throat:150mm
Executive Standard		BSEN 6506, ISO 6506, ASTM E10, GB/T231
Loading Method		Automatic (Loading/Dwell/Unloading)
Test Force Dwelling Time		1-99 S (1 second /step)
Hardness Resolution		0.1 HBW
Hardness Test Range		8-650HBW
Power Supply		220V /2A /50Hz
Dimensions(L*W*H)		Machine: 520*250*950mm Package: 720*530*1150mm
Weight		Net Weight: 164Kg Gross Weight:195kg

Standard Delivery

Name	Qty	Photo
Machine Mainframe	1set	
Φ150mm Test Anvil	1pc	
Φ80mm Test Anvil	1pc	
Φ80mm V-shape Test Anvil	1pc	
Φ10mm Hard Alloy Ball Indenter	1pc	
Φ5mm Hard Alloy Ball Indenter	1pc	
Φ2.5mm Hard Alloy Ball Indenter	1pc	
200±50HBW Brinell Hardness Block HB10/3000	1pc	

Standard Delivery

Name	Qty	Photo
100±25HBWBrinell Hardness Block HB10/1000	1set	
BRM-2Reading Microscope	1pc	
Power Cord	1pc	
Accessory Case	1pc	
Anti-dust Cover	1pc	
Instruction Manual	1copy	
Warranty Card	1copy	
Product Certificate	1copy	